


news

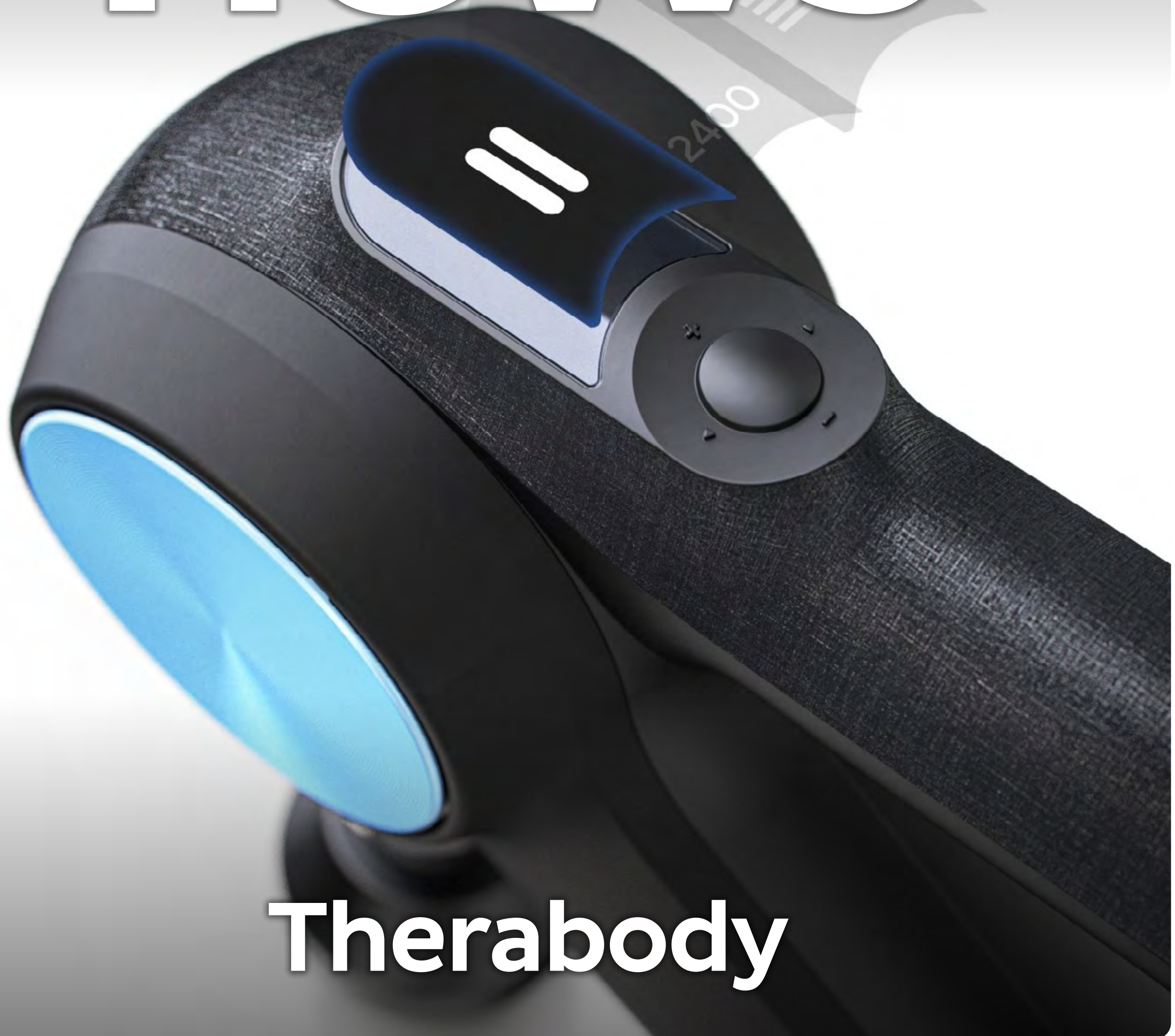

Therabody

Therabody Theragun Pro

PN **PRO-PKG-EU**
RRP **566,63 €**
EAN **810036050340**

The deep muscle treatment pros trust with the durability and features they rely on. Enhance muscle recovery, release stress and tension, and soothe discomfort with the smart percussive therapy device in a league of its own.



Therabody Theragun Pro

Smart Features

- Smart App Integration via Bluetooth®, OLED Screen, Customizable Speed Range, 3 App-Guided Preset Routines Saved on Device

Speeds

- 5 built-in speeds (1750, 1900, 2100, 2200, 2400 percussions per minute) and app-controlled customizable speed range (any speed between 1750-2400 percussions per minute).

Dimensions

- 25,5 cm x 18 cm x 7,6 cm (10 in x 7,1 in x 3 in)

Weight

- 1,3 kg (2,9 lbs)

Battery Life and Charge Time

- 300 minutes total, 150 minutes per battery; 75 minutes



Using the PRO

1. Insert the lithium-ion battery fully before turning on the PRO.
2. To turn on the PRO, firmly press and hold the center button for 2 seconds on the button pad.
3. Press the plus (+) and minus (-) buttons on the button pad to increase or decrease speed.
4. Press the left (<) and right (>) arrows on the button pad to toggle between the PRO's pre-loaded treatment presets.
5. The PRO should be used with the attachment making perpendicular contact with the skin. Only the attachment of the device should make contact with the skin or area you are treating. Do not use the device at an angle that puts the percussive arm of the device in contact with the skin.
6. To turn off the PRO, firmly press and hold the center button on the button pad until the device shuts down.

Connecting the Attachments

- Align the attachment with the tip of the rotating arm and push to connect. To remove, grip the attachment, placing 1 finger on each side of the plastic connector. Pull off the attachment.



Using the Rotating Arm

- The PRO's rotating arm adjusts to 4 positions. To adjust the arm, push the rotating arm button on the underside of the PRO and move the arm to the desired position.

Smart App Integration - Overview

- The PRO is smart device enabled with Bluetooth® connectivity, designed to pair seamlessly with the Therabody app. When you download the app, you'll be guided through a step-by-step setup and shown how to pair your device. Once paired with your device, the app can remotely adjust speed, activate guided treatment presets, and make real-time recommendations based on your usage. The PRO's OLED screen supports your app experience by displaying the current speed and force to help guide your treatment.

What's Included

- Theragun PRO
- Protective Carrying Case
- 6 Attachments with Travel Pouch - Supersoft, Dampener, Standard Ball, Wedge, Thumb, Cone
- 2 Rechargeable Lithium-ion Batteries
- Lithium-ion Battery Charger



Therabody Theragun Elite White

PN ELITE-PKG-WHT-EU
RRP 377,74 €
EAN 810036050548



The powerful deep muscle treatment, in an ultra-quiet smart percussive therapy device with advanced sound insulation. Elite works deep to melt away tension and release soreness, empowering you to care for your body daily.



Therabody Theragun Elite Black

PN ELITE-PKG-BLK-EU
RRP 377,74 €
EAN 810036050425



Therabody Theragun mini White

PN **MINI-PKG-WHT-EU**
 RRP **188,85 €**
 EAN **810036053082**



The mini is your pocket-sized partner, giving you Theragun quality muscle treatment with unparalleled portability. Compact but powerful, Theragun mini is the most agile massage device that goes wherever you do.



Therabody Theragun mini Black

PN **MINI-PKG-EU**
 RRP **188,85 €**
 EAN **810036050210**



Therabody Wave Solo

PN WAVESOLO-EU
RRP 75,89 €
EAN 810036054683



This ultra-portable smart vibration therapy device delivers pinpointed pressure and vibration to focused areas for reduced tension and improved movement. Three powerful vibration frequencies and innovative wave texture allow its ergonomic shape to work deep into hard-to-reach areas, while our innovative QuietRoll Technology™ effectively mutes sound and reverberations.



Therabody Wave Duo

PN WAVEDUO-EU
RRP 94,41 €
EAN 810036055130



Therabody Pro Wireless Charger

PN **PRO-STND-PKG-EU**
RRP **81,44 €**
EAN **810036050814**



Charge your PRO in one easy step, simply sliding the device into this convenient Wireless Charging Stand without the hassle of wires or pins. And because the PRO can remain in the stand even when fully charged, it allows you to transition from charging to using your device faster than ever.

Therabody Prime Charging Stand

PN **PRIME-STND-PKG-EU**
RRP **62,93 €**
EAN **810036051026**



Therabody Multi-Device Wireless Charger

PN **MD-CHRG-PKG-EU**
RRP **81,44 €**
EAN **810036052023**

